

A1C GO HbA1c POCT Analyzer

The A1C GO is an advanced Point-of-Care (POCT) analyzer designed for fast and accurate HbA1c testing. Utilizing Boronate Affinity Chromatography, it ensures high precision, making it an ideal choice for diabetes diagnosis and monitoring.



Key Features & Benefits



Fast & Reliable Results

Delivers precise HbA1c readings in just 2 minutes.



Boronate Affinity Chromatography Technology

No interference from hemoglobin variants, ensuring high accuracy.



Small Sample Volume

Requires just 4µL of whole blood, making it ideal for quick testing.



User-Friendly & Compact Design

Lightweight, portable, and easy-to-use touchscreen interface.




Compact



Fully
Automated



Cartridge
based

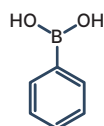

Boronate
Affinity



Add sample



Insert into Analyzer
Start testing



Technical Specifications

Measurement Technology	Boronate Affinity Chromatography
Sample Type	Whole blood (fingerstick or venous)
Sample Volume	4 μ L only
Portability	Lightweight, easy-to-use design
Test Time	2 minutes
Storage	Room temperature stable reagents
Data Connectivity	LIS/HIS system compatible

Why Choose A1C GO?

- ✓ Fast, 2-minute testing for immediate results
- ✓ High precision with no hemoglobin variant interference
- ✓ Compact, portable, and easy-to-use design